ORALSAN®

st8

Sonic
Toothbrush

Pag 2

Pag 13

Pag 24

www.ids dental.it

Dear customer, thank you for purchasing this product. Please read the user manual carefully and keep it for future use.

Product introduction

The Sonic Toothbrush adopts innovative brushless motor to produce a high vibration of 38,000 per minute. Equipped with a DuPont diamond brush head, the toothbrush possesses a superior cleaning ability and effectively removes dental plaques, food residues, stains, etc. Humanized designs that includes 5 brushing modes, auto-memory function and smart timer meet different teeth brushing needs. There is a built-in 1800mAh battery which lasts up to 50 days on a single charge. And the IPX7 waterproof technology enables you to rinse the toothbrush directly under tap. Besides, the compact appearance and ergonomically designed body also deliver a nice using experience. Continuous use of the toothbrush will bring you healthier teeth, fresher breath and more confident smile.

Important

Read this user manual carefully before you use the appliance and save it for future reference.


Warnings

1. Check if the voltage indicated on the appliance corresponds to the local mains voltage before you connect the appliance.
2. The battery may be exhausted during long-distance transportation.
A full charge will be recommended before use.
3. If the appliance is damaged in any way (brush head, brush handle, and / or charging cable), stop using it. Always have them replaced with original type in order to avoid hazard.
4. This device should not be used by persons (including children) with reduced physical, sensory or mental abilities, persons who are unfamiliar with or have no experience with the use of these devices, unless the use is under the supervision, or instructions regarding its use have been given by persons responsible for their safety.
5. It is necessary to ensure that children do not play with the device (the latter is not suitable for use by children).


Caution

1. If you have had oral or gum surgery in the previous 2 months, consult your dentist before you use the toothbrush.
2. Consult your dentist excessive bleeding occurs after using this toothbrush or if bleeding continues to occurs after 1 week of use.
3. Consult your physician if you have other medical concerns.
4. This product complies wah the mfay standards for electromagnetic devices. If you have a pacemaker or other implanted device, contact your physician or the device manufacturer prior to use.
5. This appliance has only been designed for cleaning teeth, guns md tongue. Do not use it for any other purpose. Stop using the appliance and consult your doctor if you experience any discomfort or pain.
6. The sonic toothbrush is a personal care device and is not intended for use on multiple patients in a dental practice or institution.
7. Stop using a brush head with crushed or bent bristles. Replace the brush head every 3 months or sooner if signs of wear appear. Do not use other brush heads than the ones produced or recommend by the manufacturer to avoid bad using effect or damaging the brush body.
8. If your toothpawe contains peroxide, baking soda or other bicarbonate (common in whitening toothpastes), thoroughly clean the brush head and the handle with soap and water after each use. Otherwise the plastic may crack.
9. Do not clean brush heads and the handle in the dishwasher.


Product overview

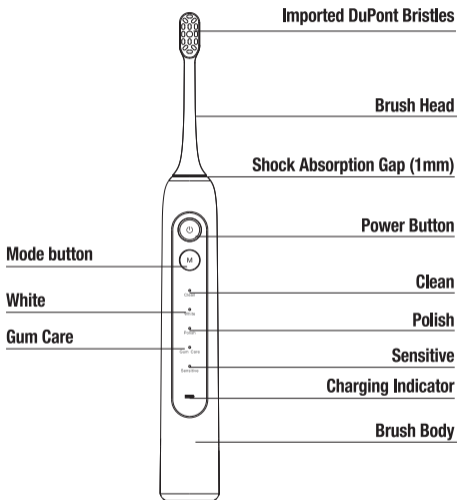
Charging method

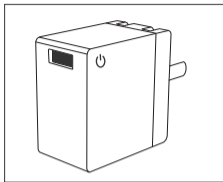
1. Prepare a qualified 5V/1A adaptor (5V/1A adaptors should be prepared by users. Most smartphone adaptors are suitable. Please make sure they are qualified before use)

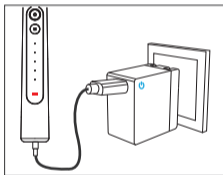2. Insert the DC port into the toothbrush and the USB port into the adaptor Then connect the power supply. The indicator will light on red .

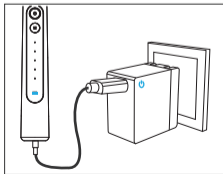
3. The indicator turns blue when fully charged, then unplug the adaptor. (A full charge takes about 6 hours.)

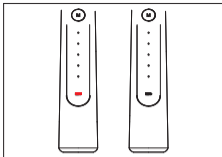
4. Rapid red flashing indicates that the battery is low, which reminds you to charge. Charging in time to ensure everyday use.


Application methods

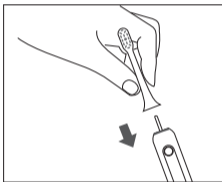
1. Install brush head: keep the direction of bristles in accordance with that of the Power Button. Press the brush head down to the metal shaft until it stops. After firmly installed, there will be a 1mm or so shock absorption gap between the brush head and the brush body (The gap is designed to make the brush head shock better and reduce the shock sense of the brush body when working, but not an issue of quality. Please feel free to use).

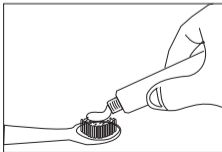
2. Squeeze toothpaste: wet the brush head and squeeze appropriate toothpaste on it.

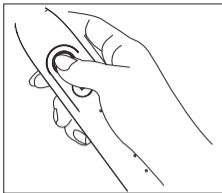
3. Press the Power Button to switch on the sonic toothbrush. Press the Mode Button to choose a mode as you prefer (The toothbrush will automatically select your last choice next time)

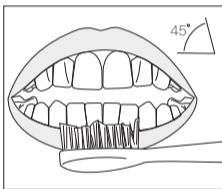
4. Put the brush head: press the brush head against the teeth at a slight angle toward the gum line.

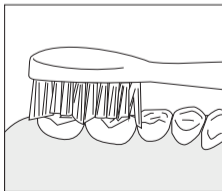
5. Apply light pressure to maximize the sonic toothbrush's effectiveness and let it do the brushing for you. The pressure should depend on your using experience.

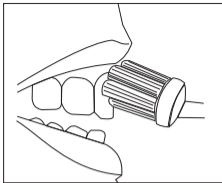6. Gently move the brush head slowly across the teeth in a small back and forth motion so the longer bristles reach between your teeth. Continue this motion throughout the brushing cycle. Note: To make sure you brush evenly throughout the mouth, divide your mouth into 4 sections.

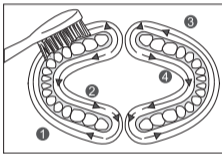7. Begin brushing in section 1 (outside upper teeth) and brush for 30 seconds before you move to section 2 (inside upper teeth). Continue brushing in section 3 (outside lower teeth) and brush for 30 seconds before you move to section 4 (inside lower teeth).

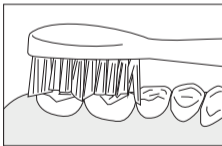
8. After you have completed the brushing cycle, you can spend additional time brushing the chewing surfaces of your teeth and areas where staining occurs.

Note: the sonic toothbrush is safe to use on:

- Braces (brush heads wear out sooner when used on braces)
- Dental restorations (fillings, crowns, veneers), choose sensitive mode

Note: you can switch off the toothbrush at any time by pressing the Power Button.

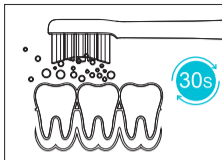9. Smart timer: each mode has a corresponding timing; remind to change brushing area every 30 seconds

10. Clean toothbrush: after brushing, please turn off the toothbrush before take it out of your mouth. Gargle and wash brush head with clean water, then shake off the water and dry brush body with a soft towel. Take care of the toothbrush.

Modes instruction

According to different teeth conditions, we designed 5 different brushing modes: Clean, White, Polish, Gum Care and Sensitive. You can choose different modes depending on your need.

1. Clean: (Standard mode) Effectively remove all kinds of dental plaque. Suitable for everyday clean of common people

2. White: 2-minute brushing effectively remove teeth stains. Extra 30-second helps polish and whiten the front teeth.

3. Polish: With intensive, high-speed and high-frequency strokes, this mode is especially designed for brightening teeth and brushing the front teeth.

4. Gum Care: Clean mode plus 1-minute to massage the gum. Helps improve gum health.


5. Sensitive: The soft vibration of the brush head gently cleans the dirt and plaque hidden in the teeth, keeps the teeth healthy and prevents gums from bleeding. Suitable for primary users, people with sensitive teeth, the elderly and adolescents.

Clean & Care

Please do not clean the brush head or the brush handle in the dishwasher, as this may cause damage.

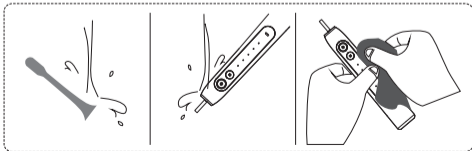1. Brush handle a. remove the brush head and wash the metal shaft area with warm water. Do not push on the rubber seal on the metal shaft with sharp object, as this may cause damage b. Wipe the whole surface of handle with a wet cloth.
2. Brush head: a. wash the brush head and bristles after each use b. remove the brush head and wash the brush head connection with warm water at least once a week.

Faq

1. Why doesn't the sonic electric toothbrush work?
It may be caused by an exhausted battery. You need to charge the toothbrush for at least 6 hours.
2. Why do I feel itching when I use the sonic electric toothbrush?
You may experience a slight itching/tingling sensation when you use the sonic toothbrush for the first time. As you become accustomed to brushing with the sonic toothbrush, this sensation diminishes.
3. The brush head chatters against my teeth. What can I do about this?
To avoid this, keep the brush head properly oriented, with the bristles placed on the gumline at a slight angle.
4. What are the advantages of a sonic electric toothbrush compared to a traditional rotating electric toothbrush?
Compared to a traditional rotating electric toothbrush, a sonic electric toothbrush makes much lower sound, causes much less abrasion to teeth and is much safer to use.
5. How should I do if I feel uncomfortable when I first use an electric toothbrush?
Some people will have such feelings when they first use an electric toothbrush. You can choose sensitive mode first and then try other modes gradually.
6. How to maintain an electric toothbrush?
There is no need for special maintenance. Please put the toothbrush in dry places after cleaning. Avoid direct sunlight and do not dry with a hair dryer.

Specifications

Product Function: Remove dental plaques and whiten teeth

Battery Capacity: 1800mAh

Charging Mode DC charging

Charging Current: DC 5V/300mAh

Charging Time: About 6 hours

Working Time: About 50 days (may vary based on different using frequency)

Bristles Material: DuPont diamond bristles

Brush Head Using Time Replace the brush head every 3 months for your health. (ifs time to replace the brush head when the blue bristles turn pale in color)

Bushing Mode Clean, White, Polish, Gum Care and Sensitive

Accessories

2 * Standard Brush Head

1 * Charging Cable

1 * User Manual

Warranty

The warranty period of the sonic toothbrush is 12 months from the date of purchase (see invoice date).

All the equipment distributed by us that will need technical interventions during or outside the warranty period, must be sent directly to the company ELCAT S.N.C who will manage the repair and invoice.

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