# ORALSAN<sup>®</sup>

## Wf7

### Desktop Water Flosser

Pag 2

Pag 15

Pag 28

[www.ids dental.it](http://www.ids dental.it)

## CONGRATULATIONS!

Thank you for selecting the Oralsan WF7 oral and nasal irrigator.

The Oralsan WF7 oral irrigator cleanses your teeth and gums, removing food particles between the teeth and reducing the bacteria that can cause gingivitis and other periodontal pathologies.

The Oralsan WF7 oral irrigator effectively cleanses and clears your sinuses gently and pleasantly, removing foreign matters, crust and other undesirable materials, and massages the cilia and nose and sinuses to their normal healthy state.

To maximize your benefit from the Oralsan WF7 oral and nasal irrigator, it is important that you read and understand the information in this instruction manual.

## CONTENTS

<b>Safety</b>	<b>3</b>
-----	
<b>Features</b>	<b>4-5</b>
-----	
<b>Operation</b>	<b>6-11</b>
-----	
<b>Maintenance</b>	<b>12</b>
-----	
<b>Troubleshooting</b>	<b>13</b>
-----	
<b>Warranty</b>	<b>14</b>
-----	

## SAFETY

When using electrical appliances, especially when children are present, precautions should be carefully followed, including, but not limited to the following:

### **READ ALL INSTRUCTIONS BEFORE USING.**

### **DANGER:**

#### **To reduce the risk of electrocution:**

- Always unplug this product immediately after use
- Do not handle charger with wet hands
- Do not place in or drop into water or other liquid
- Periodically check the cord for damage
- Do not use the unit while bathing
- Do not place or store the unit where it can fall or be pulled into a tub, sink or toilet.
- Do not reach for a unit that has fallen into water.  
Unplug immediately

## WARNING

When using electrical appliances, especially when children are present, precautions should be carefully followed, including, but not limited to the following:


#### **To reduce the of burns, electrocution, fire, or injury to persons:**

- Do not plug this device into voltage system that is different from the voltage system specified on the device.
- Never operate this product if it has a damaged, or dropped into water.
- Close supervision is necessary when this product is used by, on or near children or adults requiring assistance.
- Do not direct water under the tongue, into the ear, nose or other delicate areas.  
This appliance is capable of producing pressures that may cause serious damage in these areas.
- Use this product only for it's intended use as described in the manual or as recommended by your dental professional.



- Do not use attachments not recommended by the manufacturer.
- Fill reservoir with warm water or other mouthwash or prescribed solution as recommended by your dental professional.
- Never drop or insert objects into any opening or hose.
- Do not use while wearing any oral jewelry. Remove prior to use.
- Do not use where oxygen or aerosol sprays are being administered.
- Do not add soda, salt or any other abrasive substance to the water in the reservoir as doing so will damage your Oralsan WF7.
- If you have had oral or periodontal surgery in the previous twelve months, contact your dentist prior to use.
- Do not use this product for more than 8 minutes in each two-hour period.

## SAVE THESE INSTRUCTIONS

## FEATURES
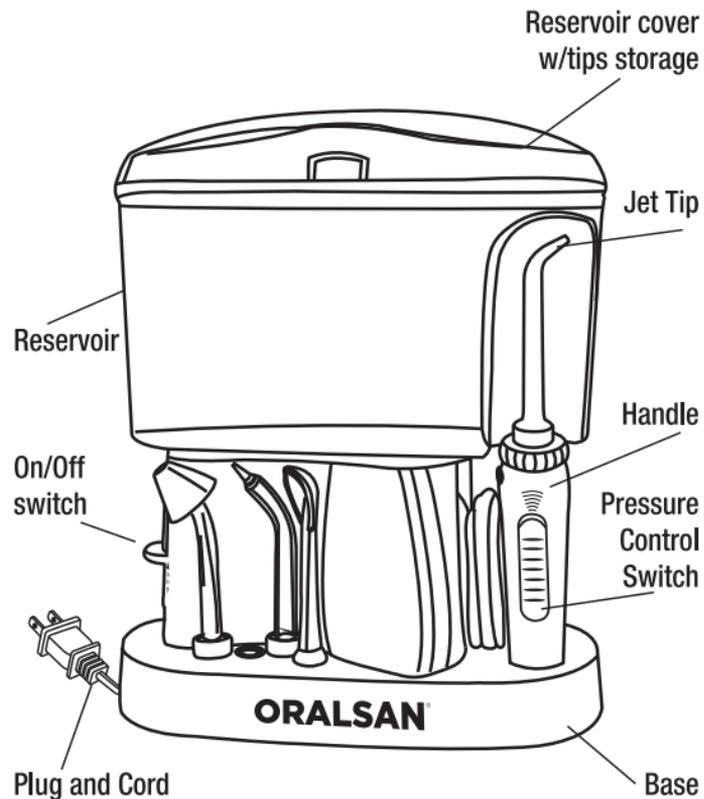
## Reservoir cover w/tips storage

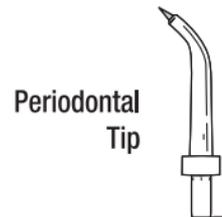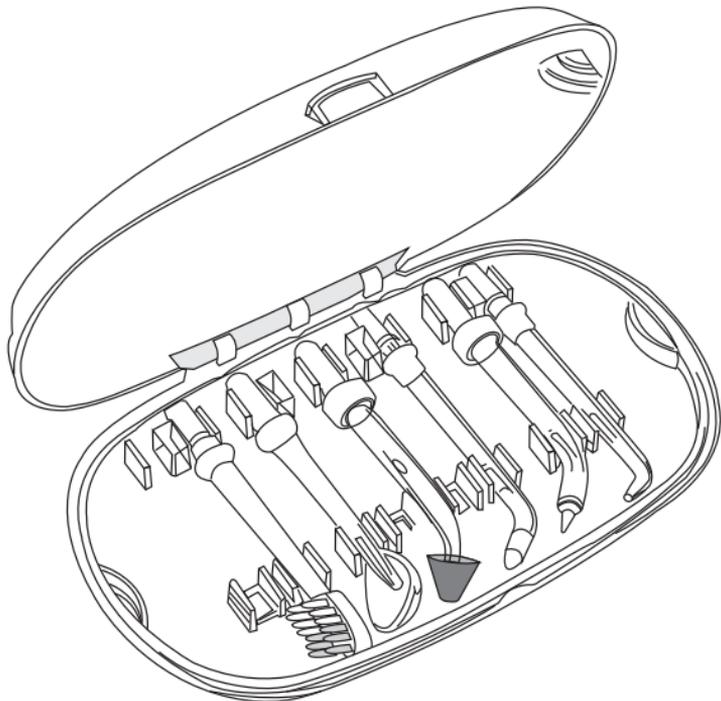
## OPERATION

Connect the device to the power supply of the voltage specified on the device.

If the power supply is connected to a wall switch, make sure it is switched on.


### Tank preparation

Lift the water tank from the base and fill with water at room temperature. Secure the tank to the base.

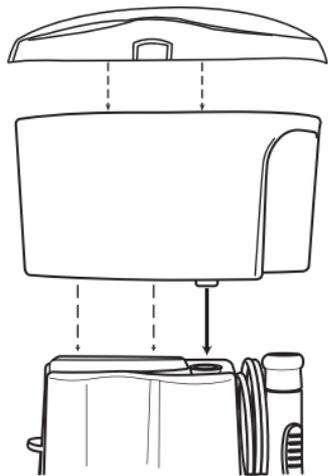
### Inserting and Removing Tips

Insert jet tip firmly into the hand. When tip is correctly inserted there should be no gap between handle and the colored tip ring. To remove tip from the handle, press the tip eject button and pull tip from the handle.
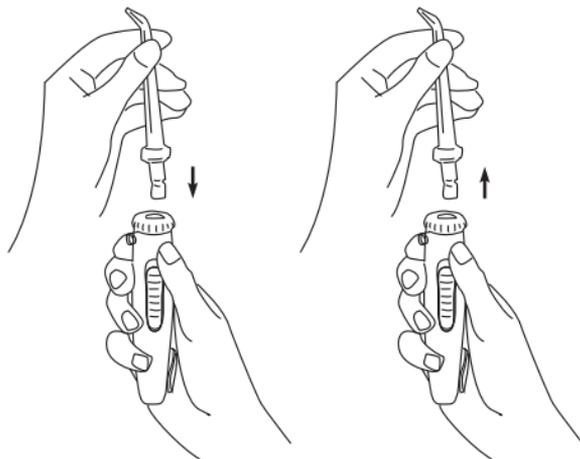
### Regulation of pressure settings

Direct the jet tip towards the sink and turn on the on / off switch in the position "On". As soon as the water begins to flow from the jet tip, move the pressure control knob to the desired setting and which is considered correct.

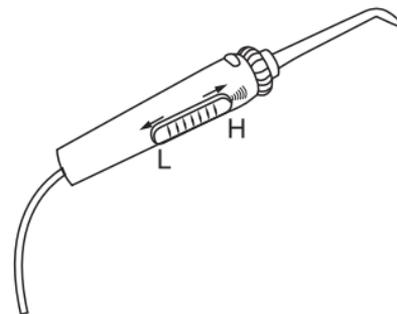
## Placement mouth


Lean low over sink and place the jet tip in mouth.

Aim the tip toward teeth.

With your free hand, slide the on/off switch on the base of the unit to the “up” position. Set pressure desired setting.

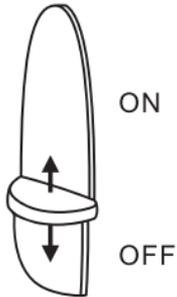Direct the stream by turning the knob at the top of the handle.
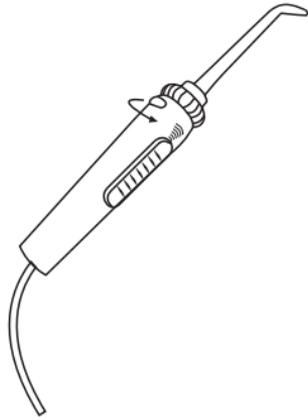
## Recommended technique


Bend over the sink and guide the jet to your teeth and gums.

Slide the handle switch up to start the water stream.

Keep your mouth slightly open to allow the water to run into the sink.

Direct the stream of water towards the gumline with at 90 degrees, so that contact is made where the tooth and gums meet.

Manoeuvre the flow of water from tooth to tooth remembering to pause briefly between each tooth.
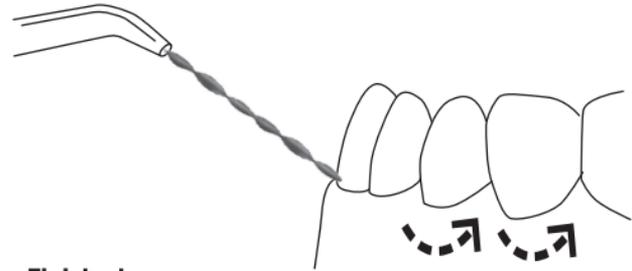
## When Finished

Turn the unit off. Empty any water remaining in the reservoir after use. Unplug the unit.

## SPECIAL TIPS


### Periodontal tip

Place the periodontal tip in contact between teeth and gumline.

Lean the soft tip against a tooth at a 45 degree angle and gently place under the gumline, into the pocket.

Start the irrigator by setting it at the lowest pressure and clean by gently following gumline, inserting the periodontal tip between teeth.

Continue tracing along the gumline to massage the gums and safely deliver solution into the gingival pockets below the gumline.
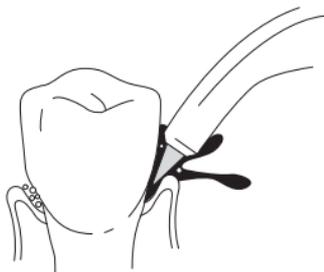
### Tongue Cleaner Tip


Lean your head over the sink. Glide the pressure button to the lowest setting. Start the irrigator and place tongue cleaner at the back of your tongue and move it forward to clean your tongue. Repeat to cover entire tongue. Set the pressure to the lowest setting and adjust to higher setting for more water pressure as it feels comfortable.

## Orthodontic tip


Use the orthodontic tip to simultaneously brush and rinse hard-to-reach areas around braces, crowns, bridges and implants.

Start on the last molar (back tooth). Place tip at 90 degree angle to the tooth at the gumline. Gently glide the tip along gumline, pausing briefly to lightly brush area between teeth and all around the orthodontic bracket, before proceeding to the next tooth and outside of both the upper and lower teeth. Continue the process until you have cleaned the inside.
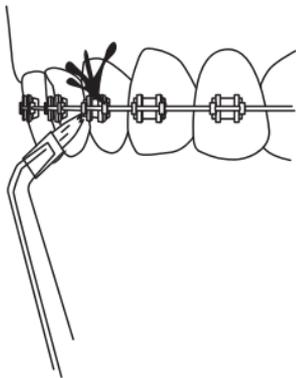
## Toothbrush Tip

Put the brush head tip into your mouth before switching on the irrigator. Move the brush head gently along your teeth and gumline region and allow the brush head to stay on the surface of each tooth for a few seconds with a slight gliding motion for the best teeth cleaning and gums massaging effect. This action loosens debris and larger food particles before flushing them away - ideal for all dental work, including fixed and removable appliances.
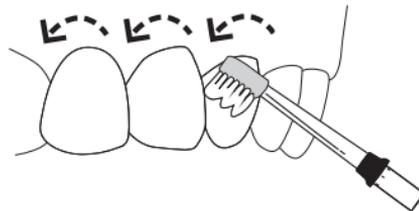
## NASAL CLEANING SYSTEM

Pulsatile nasal / sinus irrigation is an extremely effective system for cleansing the nose and sinus, an effective tool for maintenance of healthy nasal & sinus cilia.

### IMPORTANT SAVEGUARDS TO AVOID INJURY:

- Follow your physician's instructions carefully when using this product.
- Closely supervise use by children and individuals with special needs.
- Do not use the nasal tip attachment anywhere except in nose.
- Do not use any other attachment in the nose except the original nasal tip.

### WARNINGS:

- If you have an acute sinus infection, consult your doctor before using it.
- If you find yourself gagging or experience ear pain, you are trying to force too much solution into the nasal passage too quickly.
- Coughing is normal as some of the solution may run down the back of the throat.

### PREPARING RINSING SOLUTION:




For maximum comfort, use a saline solution or doctor prescribed solution.

To make the saline solution, use non-iodized or kosher, pickling, canning or sea salt. In a clean container, mix a teaspoon of salt (0.22 oz by weight) salt, ½ teaspoon (0.07 oz by weight) of baking soda (sodium bicarbonate) with 8 oz of lukewarm water. Dissolve salt thoroughly and pour into the reservoir.

### DIRECTIONS FOR USE:

- Insert the nasal jet tip into the opening on the handle. Push down firmly! Jet tip should snap into place.
- Position the soft nasal blue bulb on the sink, open end pointing upward.
- Turn "ON" the irrigator and glide the pressure setting to adjust the stream so that it is almost one inch high.

- Keep it at this pressure setting, switch “OFF” the appliance.
- Insert the soft nasal blue bulb comfortably into the nose and turn “ON” the appliance.
- Lean over the sink so you can see the drain.  
In this position, none of the solution will go in your throat, and the treatment will be comfortable and pleasant.
- The solution will gently go in one side of the nose and out the other while stimulating nasal cilia and removing undesirable material. Continue rinse until liquid that comes out clear.
- Remove nasal tip attachment from current nostril, place in other nostril and repeat the process.
- Use about 50% of the mixed solution on each side of the nose.
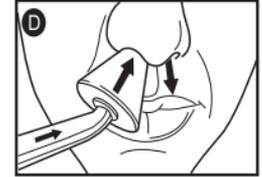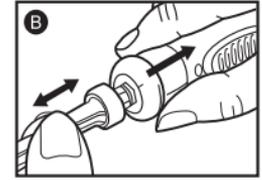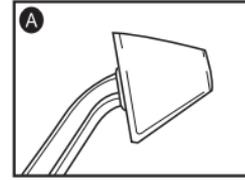
## MAINTENANCE

- Make sure to remove to throw out the water inside the tank, and run the appliance for 2 or 3 seconds to drain the water after usage. Otherwise bacteria could collect in the remaining water.
- Do not submerge the unit in water or any other liquids.
- Wipe the exterior when necessary with a soft cloth and mild non-abrasive cleanser.
- Reservoir may be washed with soapy water. Before cleaning, unplug from the electrical outlet. Wash the dirt off from inside using a brush or similar, and wipe off with a cloth.
- Dry the reservoir when not using for a long period.

## TROUBLESHOOTING

<b>PROBLEM</b>	<b>CAUSE</b>	<b>CORRECTIVE ACTION</b>
<b>Jet tips and Tongue cleaner</b> Leakage between the tip and handle knob	Tip not fully engaged	Remove tip and reinstall
<b>Periodontal tip</b> Tear in soft tip	Tip is worn	Replace with new tip
<b>Base</b> Hose to the handle leaks	Damage to the hose	Contact supplier
Reservoir leaks	Missing reservoir valve or reservoir valve is upside down	Replace reservoir valve. Turn valve right-side up
Inadequate pressure	Reservoir not positioned properly	Re-position reservoir. Press down firmly on the reservoir.
Unit does not start	Motor not running	Be sure outlet is functioning: test with another appliance (check electrical reset button and / or wall switch)

## **Accessories**

- 11 x Tips (see box)
- 1 x Cleaning Towel
- 1 x User manual

## **WARRANTY**

The supplier will repair or replace (at its option) your unit free of charge for a period of 12 months from the date of purchase if the appliance is defective in workmanship or materials.

This guarantee does not cover: damage due to improper use, normal wear or use as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorized persons and if original parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your purchase receipt to:

ELCAT S.N.C - VIA P. CRESPI, 10 - MILANO

Tel. +39 02 26111520 – fax +39 02 26143958 - e-mail: [info@elcatmilano.com](mailto:info@elcatmilano.com)

Purchaser is responsible for all postage and freight charges and any other expenses incidental to claiming against this warranty.

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