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**RECOMMENDED BY DENTISTS
AND HYGIENISTS**

For further info
ask your dentist or hygienist
or visit our website www.oralprox.it

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
9 sizes available
Ask your pharmacy for OralProx

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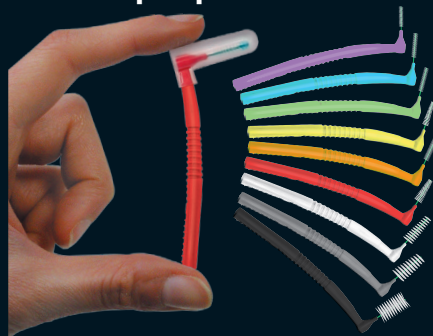**The handy pocket interdental brush
that's always ready to use!**

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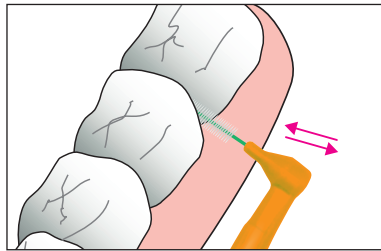
ORALPROX®

EN Why use an interdental brush

Over 80% of teeth and gum diseases occur in interdental spaces.

Food particles and bacterial plaque settle in these areas which are difficult to reach with a standard toothbrush.

Over time, these bacteria will inevitably lead to cavities and gum disease.

That's why dentists and hygienists recommend using interdental brushes.

Where to use OralProx


It is not always possible to insert an interdental brush between all teeth.

The first areas to check with OralProx are those where food tends to get trapped.

How to use OralProx

Place the brush head just above the gum line. Push the head between the teeth without using excessive force. Wash and replace the cap after use.

Using OralProx once a day is enough. You can use it more if you like.

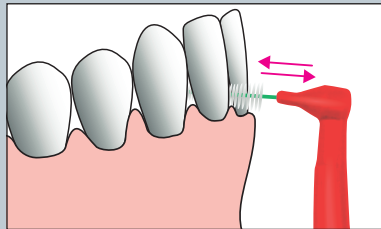ORALPROX®

What you should know

- **Bleeding:** if gums are inflamed, bleeding may occur while using the interdental brush. This should gradually reduce as the gums get better.
- **Warning:** if the brush head should break between the teeth:
 1. lean forward to avoid swallowing the fragment
 2. push it out using another OralProx
- **When to replace the brush:** to ensure that this device is always efficient and hygienic, replace OralProx at least once a week. Replace it immediately if the steel wire should become deformed.