**Waterpik® Ultra  
Cordless Dental Water Jet**  
**Waterpik® Ultra Handy Jet**  
**Hydropulseur Ultra**  
**Waterpik® sans fil**  
**Waterpik® Snoerloze**  
**Monddouche**  
**Idropulsore senza fili**  
**Waterpik® Ultra**  
**Irrigador dental inalámbrico**  
**de Waterpik® Ultra**

Congratulations! By selecting the Waterpik® Ultra Cordless Dental Water Jet you are now on your way to better gum health, and a healthier smile!


With a clinically proven combination of water pressure and pulsations, the Waterpik® Dental Water Jet removes harmful bacteria deep between teeth and below the gumline – where brushing alone won't reach.

Thank you for choosing the Waterpik® Ultra Cordless Dental Water Jet.

### Traditional Oral Care


Traditional brushing can't easily reach hidden food debris and harmful bacteria.

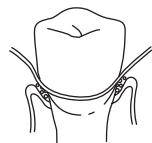Flossing has limited access below the gumline.

### Oral Care with the Waterpik® Dental Water Jet

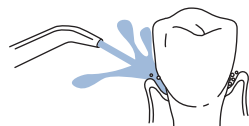The Waterpik® Dental Water Jet removes trapped food debris and harmful bacteria.

Important Safeguards	4
Product Description	6
Getting Started	8
Care of Your Waterpik® Ultra Cordless Dental Water Jet	13
Troubleshooting Guide	14
Limited Two-Year Warranty	15

## IMPORTANT SAFEGUARDS

### IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, the basic safety precautions below should always be followed.

#### READ ALL INSTRUCTIONS BEFORE USING.

#### DANGER:

##### To reduce the risk of electrocution:

- Do not handle charger with wet hands.
- Do not immerse in water or other liquid.
- Check the charger cord for damage before the first use and during the life of the product.
- Do not use while bathing.
- Do not place or store product where it can fall or be pulled into a tub, shower stall or sink.
- Do not reach for a product that has fallen into water. Unplug immediately.







### WARNING:

#### To reduce the risk of burns, electric shock, fire, serious injury or damage to tissue:

- Do not use this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into any liquid. Contact the Water Pik, Inc. service centre.
- Closely instruct and supervise children and individuals with special needs in the proper use of this product.
- Do not direct water under the tongue, into the ear, nose or other delicate areas. Without appropriate attachments recommended by Water Pik, Inc., this product is capable of producing pressures that may cause serious damage in these areas. See operating instructions for correct usage.
- Use this product only as indicated in these instructions or as recommended by your dental professional.
- Only use jet tips that are recommended by Water Pik, Inc.

- This product is to be charged by a supply which fulfills the requirements of EN61558 or equivalent Output rating as follows: 3VDC, 110mA max.
- Fill reservoir with water or other dental professional recommended solution only (i.e. mouthwash or prescribed solutions).
- Do not drop or insert any foreign object into any opening on the product.
- Keep charging unit and cord away from heated surfaces.
- Do not operate where oxygen or aerosol sprays are being used.
- Do not use mouthwashes/rinses that contain IODINE or Bleach.
- Do not use while wearing any oral jewelry. Remove prior to use.
- Do not use if you have an open wound on your tongue.
- If your physician advised you to receive antibiotic premedication before dental procedures, you should consult your dentist before using this instrument or any other oral hygiene aid.
- For household use only. Do not use outdoors.

#### SAVE THESE INSTRUCTIONS.

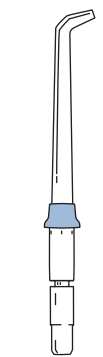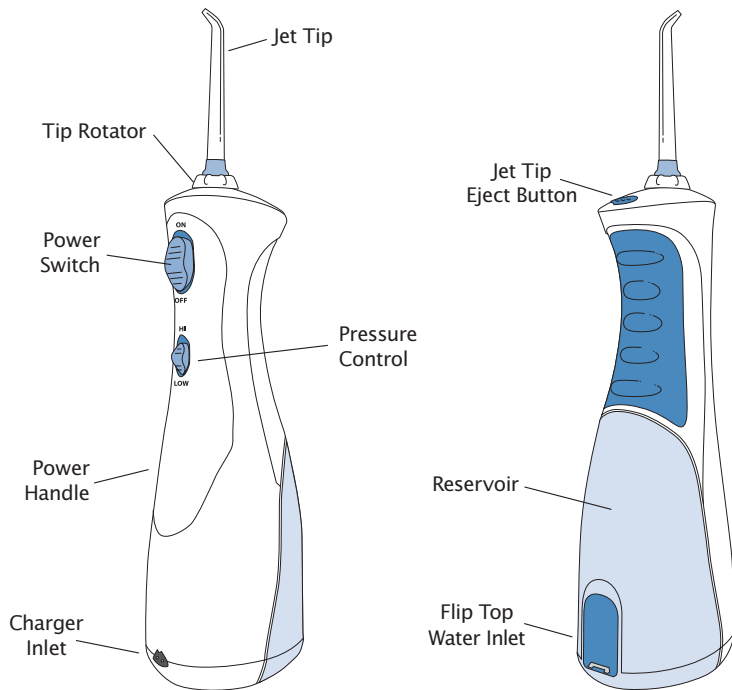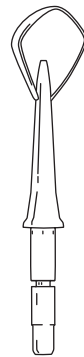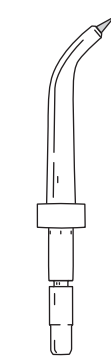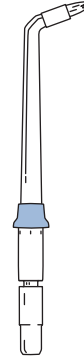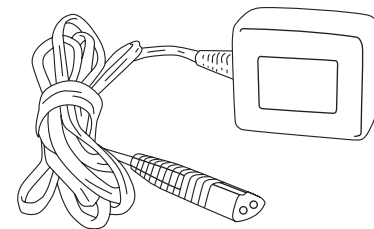Standard  
Jet Tip (1)Tongue  
Cleaner Tip (1)Pik Pocket™  
Tip (1)Orthodontic  
Tip (1)


Charger

Replacement jet tips/attachments may be purchased directly from Water Pik, Inc. call Tel +44 (0)870 734 8450.


## GETTING STARTED

**Important! Turn Dental Water Jet off and charge for 16 hours before first use.**


Plug the cord into any standard electrical outlet. The Waterpik® Ultra Cordless Dental Water Jet will stay charged for approximately 1 week depending on use. Make sure you push the switch into the 'OFF' position. The unit will not charge if the switch is in the ON position. See "Care of your Waterpik® Ultra Cordless Dental Water Jet" for guidance on how to look after your battery.

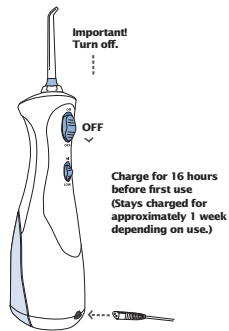**Filling the Reservoir**

Lift flip top on reservoir and fill the reservoir with lukewarm water.

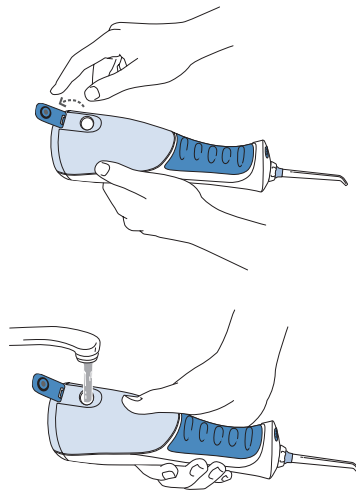**Inserting and Removing Tips**


Insert tip into the center of the knob at the top of the Dental Water Jet handle. The colored ring will be flush with the end of the knob if the tip is correctly locked into place.

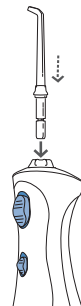
To remove the tip from the handle, press the jet tip eject button and pull tip from the handle.

**Adjusting the Pressure Setting**


Slide the pressure switch to low for first time use and for gentle cleaning of sensitive gums.

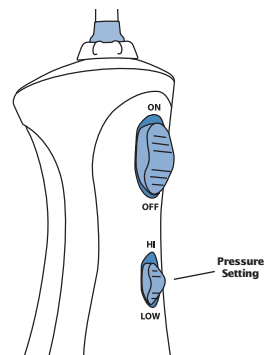Keep unit upright during use for best results.

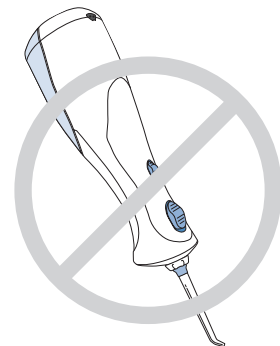
Direct the stream by turning the knob at the top of the handle.

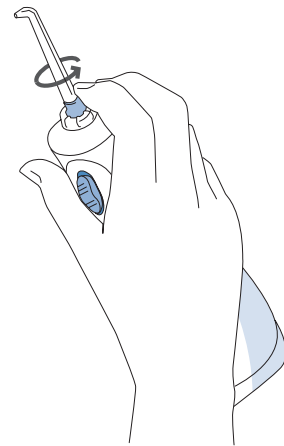
**Recommended Technique**


Lean low over sink and place the jet tip in mouth. Aim the tip toward teeth, and turn the unit on.

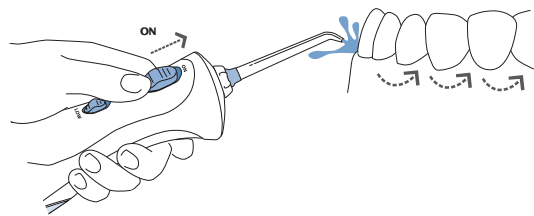
Direct the jet stream at a 90-degree angle to your gumline. **Slightly close lips to avoid splashing** but allow water to flow freely from mouth into the sink. For best results, start in the molar area (back teeth) working toward the front teeth. Follow the gumline and pause briefly between teeth. Continue until all areas around and between the teeth have been cleaned.

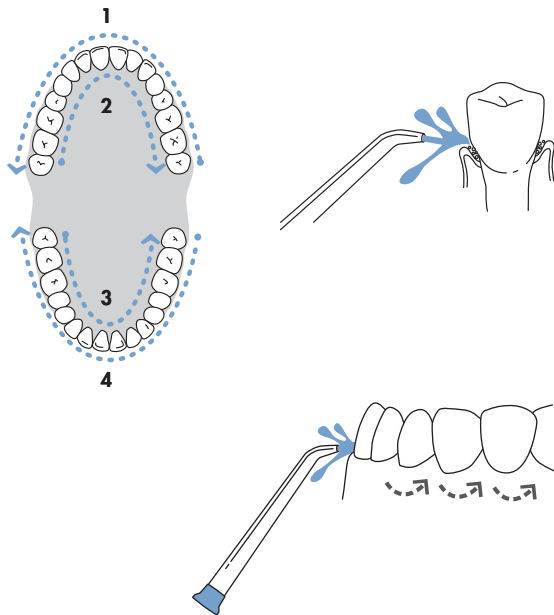**When Finished**

Turn the unit off. Empty any liquid left in the reservoir.

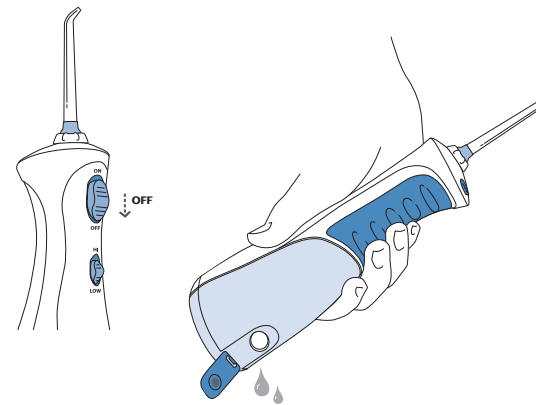If desired, the reservoir may be removed from the power handle by pulling the reservoir release latch while sliding the reservoir down towards the base of the unit.


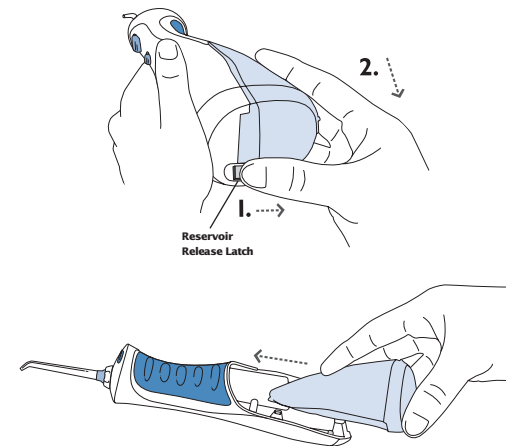
To replace the reservoir on the power handle, simply slide it up towards the top of the unit. The reservoir will snap into place when fully engaged.

## GETTING STARTED

**SPECIAL TIPS****Pik Pocket™ Tip**

The Pik Pocket™ Tip delivers water or anti-bacterial solutions deep into periodontal pockets.


To use the Pik Pocket™ Tip, lean low over your sink and place the soft tip against a tooth at a 45-degree angle and gently place the tip under the gumline, into the pocket. Turn the unit on and clean by gently following gumline, inserting tip in between teeth. Continue tracing along the gumline until the reservoir is empty.

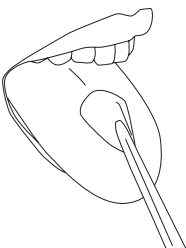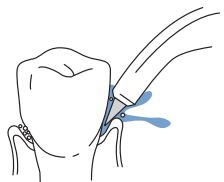**Tongue Cleaner Tip**

Place the tongue cleaner on your tongue and start the Dental Water Jet – use the lowest pressure setting, increase as comfortable. Gently place the tongue cleaner in the center/midline of your tongue about half way back and with light pressure pull forward.

**Orthodontic Tip**

Start on the last molar (back tooth). Aim tip at a 90 degree angle to the tooth at the gumline. Gently glide tip along gumline, pausing briefly to lightly brush area between teeth and all around orthodontic bracket, before proceeding to the next tooth. Continue process until you have cleaned and irrigated the inside and outside of both the upper and lower teeth.

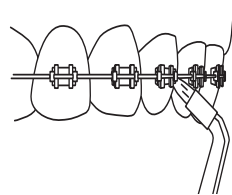For optimal results it is recommended that orthodontic tips are replaced after 6 months of use.

Contact the service center for replacement jet tips.

## CARE OF YOUR WATERPIK® ULTRA CORDLESS DENTAL WATER JET

Never immerse the unit in water. Wipe the exterior when necessary with a soft cloth and mild non-abrasive cleanser. Reservoir may be washed with soapy water or in top rack of dishwasher.

Before cleaning, unplug from the electrical outlet.

Before exposing the Dental Water Jet to freezing temperatures, remove the reservoir and run the unit until completely empty.

To maximize the capacity and life of the rechargeable batteries, the unit should be charged for 16 hours before the first use. Then unplug the charger and do not recharge until the unit has been fully discharged by daily use. The charger may then be reconnected and left plugged in when not in use. For best results, the unit should be completely discharged by use twice a year, and then fully recharged for 16 hours after each discharge. Make sure you push the switch into the 'OFF' position. The unit will not charge if the switch is in the ON position.

**Battery disposal**

The unit contains a Nickel Cadmium battery. At the end of the battery life, please recycle the unit at your nearest recycle center. The battery is not replaceable. For unit availability, contact your nearest Authorised Service Center.

**Service Maintenance**

Waterpik® Dental Water Jets have no consumer-serviceable electrical items and do not require routine service maintenance. Repair of this product should be performed only by an authorized service center or by Water Pik, Inc. For all your service needs go to [www.waterpik.com](http://www.waterpik.com).

Some replacement parts and accessories for the Ultra Cordless are also available for sale at the customer service number shown below.

Refer to the serial and model numbers in all correspondence. These numbers are located on the bottom of your unit.

UK Service Centre, call Tel +44 (0)0870 734 8450.

## TROUBLESHOOTING GUIDE

PROBLEM	CAUSE	CORRECTIVE ACTION
<b>Jet Tips</b> Leakage between the tip and the tip rotator.	1) Tip not fully engaged.	1) Remove tip and reinstall.
<b>Dental System</b> Reservoir leaks.	1) Liquid was left in reservoir. 2) Base not correctly engaged with reservoir. 3) Water inlet is open.	1) Always empty the reservoir after use. 2) Remove reservoir then reinstall as per instructions. 3) Close water inlet.
Inadequate pressure/No fluid flow.	1) Unit is upside down. 2) Insufficient liquid in reservoir. 3) Reservoir connection clogged. 4) Pressure setting on low.	1) Jet tip must be higher than reservoir for correct water flow. 2) Fill with liquid. 3) Fill reservoir with warm tap water and one cap of vinegar. Run until empty. Rinse with full reservoir of tap water. 4) Adjust pressure setting to high.
Unit does not charge	1) Mains power socket is not live 2) Switch is in ON position during charging	1) Turn ON the mains power socket; try an alternative mains power socket. 2) Turn unit switch to OFF position for charging.
Unit doesn't start.	1) Battery not charged.	1) Be sure the power socket is functioning; test with another appliance. Turn unit off. Plug in charger. Charge for 16 hours. 2) Be sure charger pins are clean and unobstructed.

Still have questions? Go to [www.waterpik.com](http://www.waterpik.com)

If none of these corrective actions are successful, call the Water Pik, Inc. Service Centre +44 (0)0870 734 8450.

## LIMITED TWO-YEAR WARRANTY

Water Pik, Inc. warrants to the original consumer of this new product that it is free from defects in materials and workmanship for 2 years from the date of purchase. Consumer will be required to submit the original purchase receipt as proof of purchase date and if requested, the entire product, to support a warranty claim. Water Pik, Inc. will repair or replace at its discretion any part of the product, which in its opinion is defective, provided the product has not been abused, misused, altered or damaged after purchase, was used according to instructions, and was used only with accessories or consumable parts approved by Water Pik, Inc. Installation is the responsibility of the consumer and is not covered by the warranty. This limited warranty excludes accessories or consumable parts such as jet tips, etc.

This warranty gives you specific legal rights. You may also have other rights that vary, depending on the law in your place of residence.

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©2007 Water Pik, Inc.
Waterpik® and Waterpik® (stylized) are registered trademarks of Water Pik, Inc.
treat yourself better™ and Pik Pocket™ are trademarks of Water Pik, Inc.
Waterpik® und Waterpik® (stilisiert) sowie sind eingetragene Warenzeichen von Water Pik, Inc.
Waterpik® et Waterpik® (stylisé) sont des marques déposées de Water Pik, Inc.
Waterpik® en Waterpik® (gestileerd) zijn geregistreerde handelsmerken van Water Pik, Inc.