
waterpik®
treat yourself better™

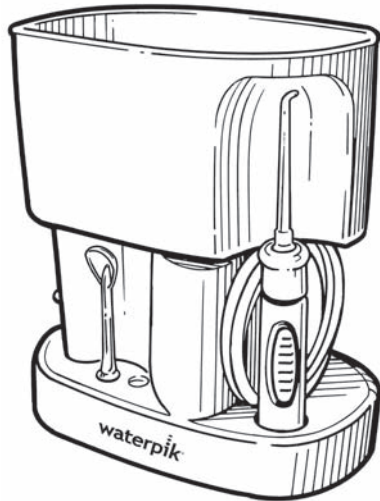Waterpik® Dental Water Jets
Waterpik® Mundduschen
Hydropulseurs Waterpik®
Idropulsori Waterpik® Dental Water Jet
Waterpik® dentale waterjets
Irrigador bucal Waterpik®

.....
www.waterpik.com

WELCOME

Congratulations! By selecting the Waterpik® Dental Water Jet you are now on your way to better gum health, and a healthier smile!

With a clinically proven combination of water pressure and pulsations, the Waterpik® Dental Water Jet removes harmful bacteria deep between teeth and below the gumline – where brushing alone won't reach.

Thank you for choosing the Waterpik® Dental Water Jet.

Traditional Oral Care


Traditional brushing can't easily reach hidden food debris and harmful bacteria.

Flossing has limited access below the gumline.

Oral Care With The Waterpik® Dental Water Jet

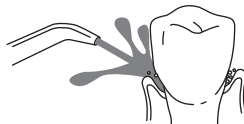The Waterpik® Dental Water Jet removes trapped food debris and harmful bacteria.

TABLE OF CONTENTS

Important Safeguards	4
Product Description	6
Getting Started	8
Care of Your Waterpik® Dental Water Jet	13
Troubleshooting Guide	14
Limited Two-Year Warranty	15

IMPORTANT SAFEGUARDS

IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, the basic safety precautions below should always be followed.

READ ALL INSTRUCTIONS BEFORE USING.

DANGER:

To reduce the risk of electrocution:

- Always unplug product after using.
- Do not handle plug with wet hands.
- Do not place or store product where it can fall or be pulled into a tub, shower stall or sink.
- Do not use while bathing.
- Do not immerse in water or other liquid.
- Do not reach for a product that has fallen into water. Unplug immediately.
- Check the power supply cord for damage before the first use and during the life of the product.

WARNING:

To reduce the risk of burns, electric shock, fire, serious injury or damage to tissue:


- Do not use this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into any liquid. Return the product to an authorized service center or Waterpik Technologies.
- Closely instruct and supervise children and individuals with special needs in the proper use of this product.
- Do not direct water under the tongue, into the ear, nose or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas. See operating instructions for correct usage.
- Use this product only as indicated in these instructions or as recommended by your dental professional.
- Do not use this product for more than 5 minutes in each two-hour period.


- Only use attachments, including jet tips, Pik Pocket® tips, tongue cleaner tips and orthodontic tips that are recommended by Waterpik Technologies.
- Do not remove the Pik Pocket® flexible tip from the clear shaft. If it separates for any reason, discard both the tip and shaft and replace with a new Pik Pocket® Tip.
- Fill reservoir with water or other recommended solutions only.
- Do not drop or insert any foreign object into any opening or hose.
- Keep electrical cord away from heated surfaces.
- Do not operate where oxygen or aerosol sprays are being used.
- Do not use outdoors.

- Do not use while wearing any oral jewelry. Remove prior to use.
- Do not use if you have an open wound on your tongue.
- If your physician or cardiologist has advised you to receive antibiotic premedication before dental procedures, you should consult your dentist before using this instrument or any other oral hygiene aid.


SAVE THESE INSTRUCTIONS.

PRODUCT DESCRIPTION

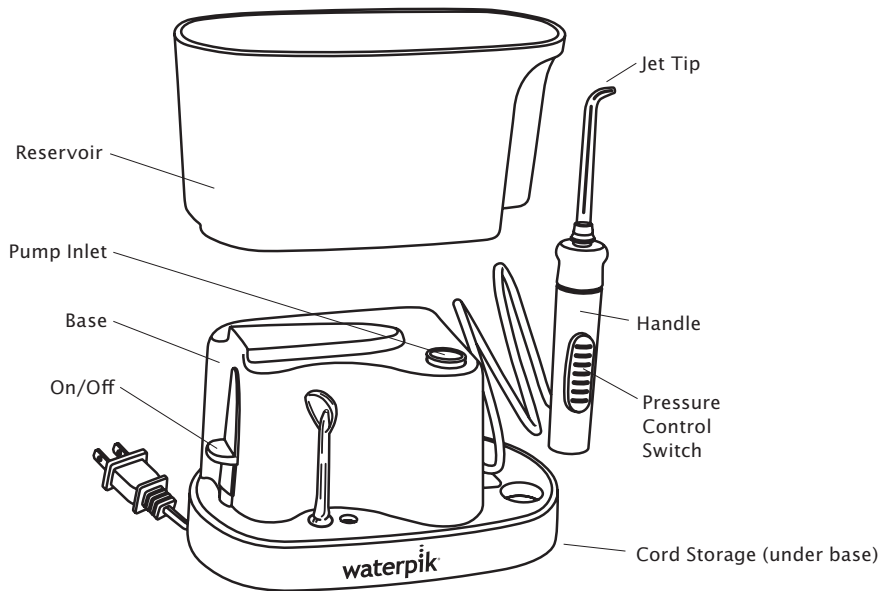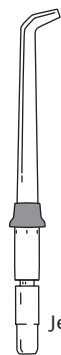
Jet Tip

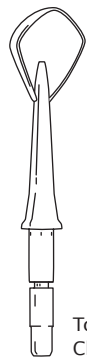Tongue
Cleaner Tip

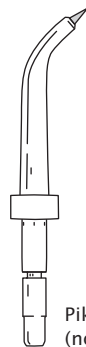
Pik Pocket® Tip
(not standard in
all models)

GETTING STARTED

Plug the cord into an electrical outlet. If the outlet is controlled by a wall switch, be sure it is turned on.


Preparing the Reservoir

Remove reservoir from base. Fill the reservoir with lukewarm water. Replace the reservoir over the base and press down firmly.

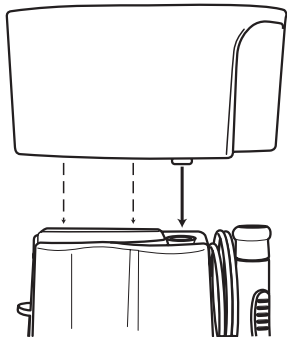
Inserting and Removing Tips

Pull the knob on the top of the handle toward the hose and insert tip into the center of the knob. The colored ring will be flush with the end of the knob if the tip is correctly locked into place. To remove tip from the handle, pull the knob on the top of the handle toward the hose and pull tip from handle.

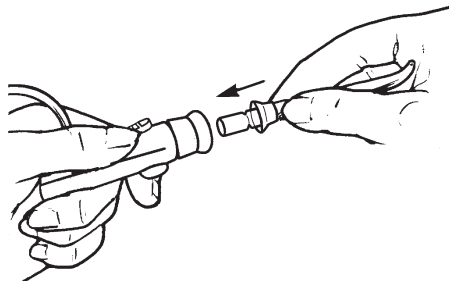Adjusting the Pressure Setting


Slide the pressure control switch located on the handle to the lowest setting (towards the tip as far as possible) for first-time use. Gradually increase pressure over time to the setting you prefer – or as instructed by your dental professional.

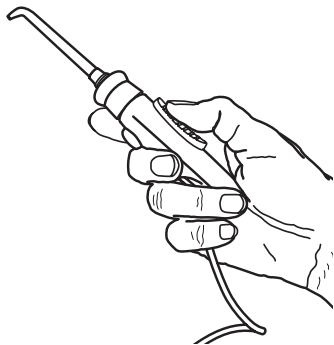
Placement in Mouth

Lean low over sink and place the jet tip in mouth. Aim the tip toward teeth.

With your free hand, slide the on/off switch on the base of the unit to the “up” position. Set pressure at desired setting.


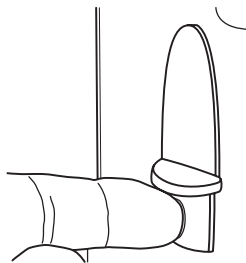
Direct the stream by turning the knob at the top of the handle.

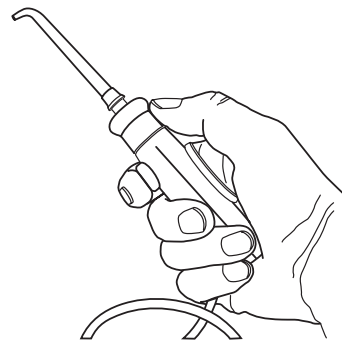GETTING STARTED


Recommended Technique

Direct the jet stream at a 90-degree angle to your gumline. Close lips to avoid splashing but allow water to flow freely from mouth into the sink. For best results, start in the molar area (back teeth) working toward the front teeth. Glide jet tip along gumline and pause briefly between teeth. Continue until all areas around and between the teeth have been cleaned.

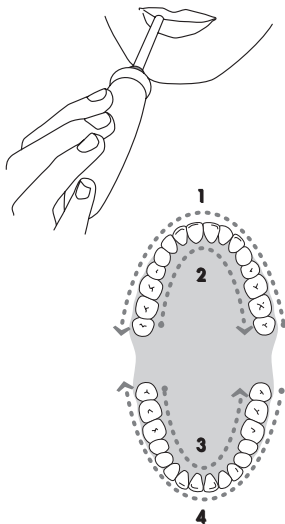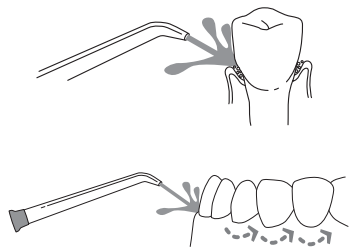Pause Control

You can temporarily stop the flow any time by pressing the pause button on the handle.

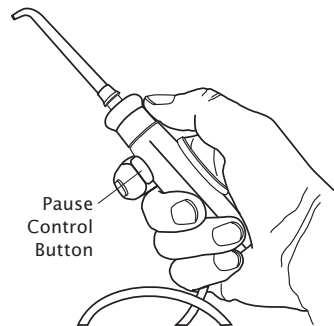When Finished



Turn the unit off. Empty any liquid left in the reservoir. Unplug the unit.

SPECIAL TIPS

Pik Pocket® Tip


The Pik Pocket® Tip delivers water or anti-bacterial solutions deep into periodontal pockets.

To use the Pik Pocket® Tip, lean low over your sink and place the soft tip against a tooth at a 45-degree angle and gently place the tip under the gumline, into the pocket. Turn the unit on to the lowest pressure setting and clean by gently following gumline, inserting tip in between teeth. Continue tracing along the gumline until the reservoir is empty.

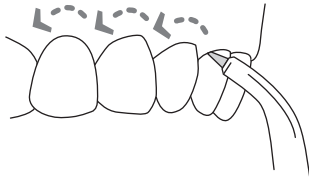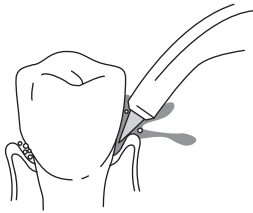Tongue Cleaner Tip

Lean over the sink. Place the tongue cleaner on your tongue and start the Dental Water Jet – use the lowest pressure setting, increase as comfortable. Gently place the tongue cleaner in the center/midline of your tongue about half way back and with light pressure pull forward.

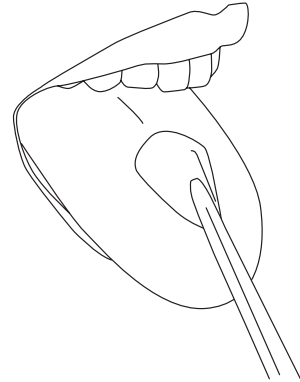GETTING STARTED

Using Mouthwash and Other Solutions

Your Waterpik® Dental Water Jet can be used to deliver mouthwash and antibacterial solutions. Ask your dental professional to recommend the appropriate solution for you. Note: use of some antibacterial solutions could shorten the life of your Waterpik® Dental Water Jet. If you have questions about using a particular solution in your Dental Water Jet, contact Waterpik Technologies customer service. After using any special solution in the reservoir, take a moment after every use to rinse the unit. Partially fill the reservoir with warm tap water. Aim the jet tip into your sink. Turn on the unit and let it run until empty. This will help keep the tip from clogging and help your unit to last longer.

CARE OF YOUR WATERPIK® DENTAL WATER JET

Clean the product when necessary by using a soft cloth and mild non-abrasive cleanser.

Before cleaning, unplug from the electrical outlet.

Before exposing the Dental Water Jet to freezing temperatures, remove the reservoir and run the unit until completely empty.

Service Maintenance

Waterpik® Dental Water Jets have no consumer-serviceable electrical items and do not require routine service maintenance. Repair of this product should be performed only by an authorized service center or by Waterpik Technologies. For all your service needs go to www.waterpik.com and click on “product support” and “service locator” for your nearest service center. Refer to the serial and model numbers in all correspondence. These numbers are located on the bottom of your unit.

When outside the United States or Canada: Contact your local dealer or write Waterpik Technologies for the address of the nearest authorized service center.

Replacement jet tips/attachments are available from select retail stores or you may purchase them directly from Waterpik Technologies.

TROUBLESHOOTING GUIDE

PROBLEM	CAUSE	CORRECTIVE ACTION
Jet Tips and Tongue Cleaner Leakage between the tip and handle knob.	Tip not fully engaged.	Remove tip and reinstall.
Pik Pocket® Tip Tear in soft tip.	Tip is worn.	Replace with new tip.
Dental Water Jet Base Hose to the handle leaks.	Damage to the hose.	Can be replaced by authorized service center or contact Waterpik Technologies Customer Service Center for a self service repair kit.
Reservoir leaks.	Missing reservoir valve or reservoir valve is upside down.	Replace reservoir valve. (Call Waterpik Technologies Customer Service Center). Turn valve right-side up.
Inadequate pressure.	Reservoir not positioned properly.	Re-position reservoir. Press down firmly on the reservoir.
Unit doesn't start.	Motor not running.	Be sure outlet is functioning; test with another appliance (check electrical reset button and/or wall switch).

LIMITED TWO-YEAR WARRANTY

Waterpik Technologies warrants to the original consumer of this new product that it is free from defects in materials and workmanship for 2 years from the date of purchase. Consumer will be required to submit the original purchase receipt as proof of purchase date and if requested, the entire product, to support a warranty claim. Waterpik Technologies will repair or replace at its discretion any part of the product, which in its opinion is defective, provided the product has not been abused, misused, altered or damaged after purchase, was used according to instructions, and was used only with accessories or consumable parts approved by Waterpik Technologies. Installation is the responsibility of the consumer and is not covered by the warranty. This limited warranty excludes accessories or consumable parts such as jet tips, etc.

This warranty gives you specific legal rights. You may also have other rights that vary, depending on the law in your place of residence.

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©2006 Water Pik, Inc.
Waterpik®, Waterpik (stylized)® and Pik Pocket® are registered trademarks of Water Pik, Inc.
treat yourself better™ is a trademark of Water Pik, Inc.
Form No. 20009170-F AB